


NURSES AND ANCILLARY STAFF

Navigating Interactive View

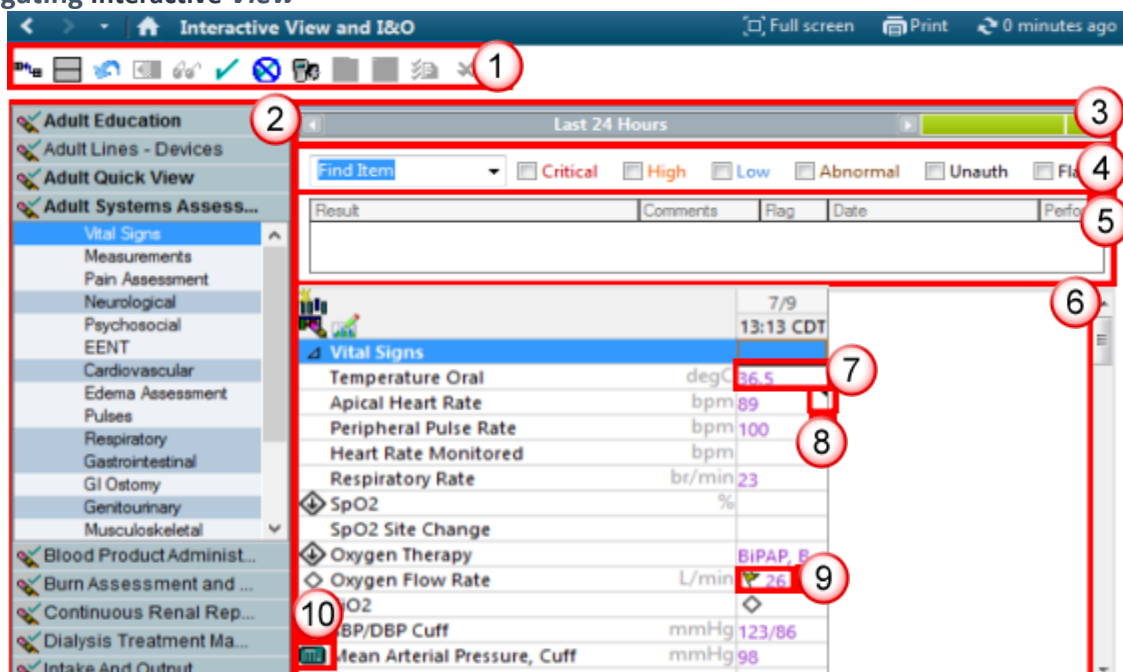

The screenshot shows the 'Interactive View and I&O' interface. The top toolbar (1) includes icons for home, full screen, print, and refresh. The left sidebar (2) is a 'Navigator' with categories like 'Adult Education', 'Vital Signs', and 'Blood Product Administ...'. The top right shows a 'Date Range Bar' (3) set to 'Last 24 Hours'. Below it are 'Filter Options' (4) for 'Critical', 'High', 'Low', 'Abnormal', 'Unauth', and 'Flag'. A 'Filter Window' (5) displays search results. The main area (6) shows a table of results for 'Vital Signs' with columns for 'Result', 'Comments', 'Flag', 'Date', and 'Perfo'. A specific result for 'Temperature Oral' (36.5) is highlighted (7). A comment icon (8) is next to it. A flag icon (9) is next to 'Oxygen Flow Rate'. A calculation icon (10) is next to 'BP/DBP Cuff'.

- 1. Toolbar:** Allows you to take action in Interactive View. See the Toolbar on page 2 for more information.
- 2. Navigator:** Displays available views. Views are displayed as bands and the sections available in each view are displayed under the bands.
- 3. Date Range Bar:** Displays date results. Use the date range bar to select the retrieval range for results.
- 4. Filter Options:** Displays filters to narrow your results. Select the appropriate filters to narrow your search.
- 5. Filter Window:** Displays search results.
- 6. Interactive View:** Displays results, the selected view, and the date in the main window in Interactive View.















- 7. Cell:** Displays results and charted values. Double-click a cell and enter a result to chart a result. Enter the result then click the **sign** button on the toolbar.
- 8. Comment Icon:** Allows you to view a result comment. Right-click a result and select **Add Comment** to add a comment.
- 9. Flag Icon:** Displays results that are flagged. Right-click a result and select **Flag** to flag a result.
- 10. Calculation Icon:** Displays calculated results. Right-click the result with a calculation icon to view the calculations.



INTERACTIVE VIEW

HOSPITAL INFORMATION SYSTEM (HIS)

Toolbar

	Collapse Navigator: Hides and displays the navigator.
	Launch Fundal Height Graph: Displays the Fundal Height graph.
	Split Window-Remove Split: Adds or removes the view to the right of the navigator, displaying multiple interactive views.
	Show Empty Columns/Rows: Displays and hides all columns and rows without results.
	Review Results: Marks all new results visible as Reviewed.
	Sign: Signs all documented results.
	Cancel: Cancels documentation without saving changes.
	Launch Labor Graph: Opens the Labor graph.
	Authenticate: Authenticates results.
	Associate Monitor: Displays the Associate Monitor window for a medical device, to disassociate, associate, or change the monitor associated with a patient.
	Add Order: Displays the order window to place a new order for the patient.
	Orders for Signature: Displays the orders scratchpad, listing the orders that need to be signed.
	Actions: Displays a list of actions you can take for a result.
	Reset Navigator Bands: Will reset the navigator bands that are displayed.

Flagging Annotations

1. Select the date and time column.
2. Click **Action** and select **Flag Annotation**.
3. Enter the date, time, title, and annotation.
4. Click **Sign**.

Uncharting a Result

1. Right-click the result you want to unchart and select **Unchart**.
2. If applicable, enter a comment.
3. Click **Sign**.

Flagging Patient Events

1. Select then right-click a result.
- Tip:** To select multiple results at once, hold **CTRL**

- as you select the appropriate results.
2. From the menu, select **Flag**.
 3. If you want to leave a comment, select **Flag with Comment** then enter a comment.
 4. Click **OK**.

Retroactively Associating a BMDI Monitor With a Patient

1. Click **Associate Monitor** and select **Modify**.
2. Select the Association Date/Time column for the patient and BMDI monitor combination you want to update.
3. Click **View Acquired Data**. Review the data from the device to determine the correct retroactive date and time to select.
4. Enter the new Association Date/Time or select a time column.

Note: If the patient you want to retroactively associate a monitor with is already associated with a BMDI monitor, you must disassociate that patient and then associate again them with the retroactive date and time selected.



INTERACTIVE VIEW


HOSPITAL INFORMATION SYSTEM (HIS)

Documenting Conditional Logic Tasks

1. Document a trigger response for a DTA with the trigger icon.
2. Document responses in the conditional fields.
3. Click **Sign**.


Adding and Defining a Dynamic Group

1. Click **Add a Repeatable Group**.

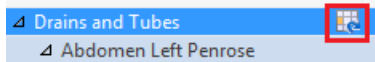

2. Select the label placeholder then select the appropriate options to define the label.
Note: The label supports up to 100 characters.
3. Click **OK**.

Inactivating a Dynamic Group Section

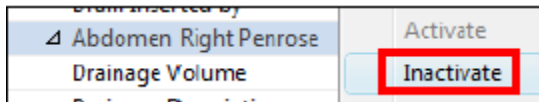

Right-click the dynamic group